



## Speed, Performance, Value - The Perfect Balance


The Adam Anti-Vibration Table helps provide a stable surface and reduces vibrations, shock and outside influences effecting devices placed on the marble working surface. Designed for the laboratory or metrology room the anti-vibration table is ergonomically styled yet built for durability.

- Large working surface of 400mm x 450mm / 16 x 18"
- Easy to set up light weight aluminium frame
- Adjustable levelling feet for uneven surfaces
- Scratch resistant marble slab
- Polished finish

The table comprises two separate parts:

- An internal working surface consisting of a large marble slab mounted on adjustable rubber shock absorbing mounts.
- An external service table with aluminium frame and adjustable levelling feet.

Overall Dimensions (w x d x h)	800 x 600 x 782 mm / 32 x 24 x 31in
Work Surface Dimensions	400 x 450 mm / 16 x 18in
Weight of Work Surface	49kg / 108Lb
Total Weight	77kg / 170Lb

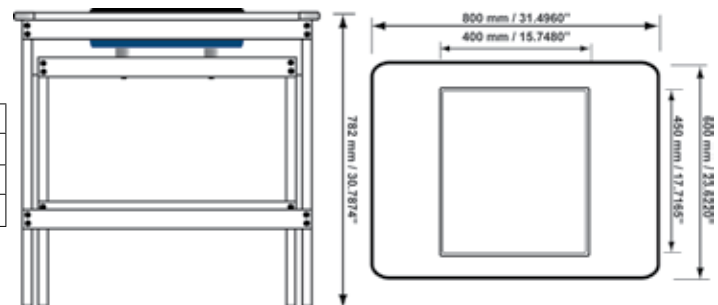

[www.adamequipment.com](http://www.adamequipment.com)

**Adam Equipment - offices worldwide**